



FLINN
SCIENTIFIC, INC.
"Your Safer Source for Science Supplies"

P.O. Box 219 • Batavia, IL 60510
(800) 452-1261 • Fax (866) 452-1436
www.flinnsci.com • E-mail: flinn@flinnsci.com
© 2012 Flinn Scientific, Inc. All Rights Reserved.

Publication No. 10216

Sudsy Kinetics

A Foamey Chemical Demonstration

Introduction

Teach kinetics concepts in a fun and sudsy way! This demonstration provides an interesting twist on the traditional "Old Foamey" or "Elephant Toothpaste" reaction. Not only will your students be amazed at the sudsy eruption—they will learn kinetics concepts along the way!

Concepts

- Kinetics/Catalysts
- Reaction intermediates
- Decomposition reactions
- Test for oxygen gas

Materials

- | | |
|--|--|
| Alconox [®] detergent, 3–4 g | Graduated cylinder, 500-mL |
| Hydrogen peroxide, H ₂ O ₂ , 30%, 10%, 3%, 20 mL of each | Large, plastic demonstration tray, several inches deep |
| Sodium iodide solution, NaI, 2 M, 4–5 mL | Lighter or matches and wood splint |
| Graduated cylinders, 10-mL and 100-mL, 3 of each | Spoon or scoop |

Safety Precautions

Hydrogen peroxide solution, 30%, is a strong oxidizing agent and a dangerous fire and explosion risk; it is severely corrosive to the skin, eyes, and respiratory tract. Do not heat this substance. Sodium iodide solution is slightly toxic by ingestion. Do not stand over the reaction; steam and oxygen are produced quickly. Wear chemical splash goggles, chemical-resistant gloves, and a chemical-resistant apron. Please review current Material Safety Data Sheets for additional safety, handling, and disposal information.

Procedure

Part 1. Effect of Concentration on the Rate of the Reaction

1. Place three 100-mL graduated cylinders in a large, plastic demonstration tray. Add 20 mL of 30% hydrogen peroxide to the first cylinder, 20 mL of 10% hydrogen peroxide to the second cylinder, and 20 mL of 3% hydrogen peroxide to the third cylinder.
2. Add 1 small scoop (3–4 g) of solid Alconox[®] detergent to each cylinder and swirl to dissolve the detergent.
3. Measure out 5 mL of 2 M sodium iodide solution in each of three 10-mL graduated cylinders. Ask your students to predict the rate at which each of the peroxide solutions will react with the iodide.
4. Ask for three student volunteers. Make sure the students are wearing chemical splash goggles; warn them to step back as soon as they pour. Have the students simultaneously pour the sodium iodide solution into the three cylinders containing the differing concentrations of hydrogen peroxide. Make observations. White foam erupts from the cylinder with the 30% peroxide the fastest, the 10% peroxide next, and only slowly rises up from the cylinder with 3% peroxide.

Part 2. Old Foamey—Observing a Reaction Intermediate and Products

1. Place a 500-mL graduated cylinder in a large, plastic demonstration tray. Measure out 20 mL of 30% hydrogen peroxide and add it to the cylinder.

© 2012, Flinn Scientific, Inc. All Rights Reserved. Reproduced for one-time use with permission from Flinn Scientific, Inc. Batavia, Illinois, U.S.A. No part of this material may be reproduced or transmitted in any form or by any means, electronic or mechanical, including, but not limited to photocopy, recording, or any information storage and retrieval system, without permission in writing from Flinn Scientific, Inc.

CHEM-FAX™ . . . makes science teaching easier.

10216
033009

2. Add 1 small scoop (3–4 g) of solid Alconox® detergent to the cylinder and swirl the mixture to dissolve the detergent.
3. Measure out 5 mL of 2 M sodium iodide solution and, quickly but carefully, pour this into the cylinder. In a few seconds, copious amounts of white foam will be produced. Observe closely at the beginning of the reaction. A brown foam is produced at first but then turns white before it erupts out of the cylinder. This is due to the presence of the free iodine produced by the extreme oxidizing ability of the 30% hydrogen peroxide.
4. Notice the steam coming off the foam—this indicates that the decomposition reaction is quite exothermic. Light a wood splint and blow out the flame. Insert the glowing wood splint into the foam. The wood splint will re-ignite in the foam—this indicates that the gas in the foam is pure oxygen. Take the glowing splint out of the foam, re-insert it, and watch it re-ignite again. This can be repeated numerous times.

Disposal

Please consult your current *Flinn Scientific Catalog/Reference Manual* for general guidelines and specific procedures governing the disposal of laboratory waste. The foam and any solution left in the cylinder or on the plastic tray may be rinsed down the drain with excess water according to Flinn Suggested Disposal Method #26b.

Connecting to the National Standards

This laboratory activity relates to the following National Science Education Standards (1996):

Unifying Concepts and Processes: Grades K–12

Evidence, models, and explanation

Content Standards: Grades 5–8


Content Standard B: Physical Science, properties and changes of properties in matter

Content Standards: Grades 9–12

Content Standard B: Physical Science, structure and properties of matter, chemical reactions, interactions of energy and matter

Discussion

Hydrogen peroxide decomposes to produce oxygen and water according to the decomposition reaction shown below:

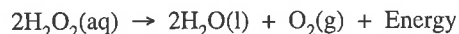

The reaction is quite slow unless catalyzed by a substance such as iodide ions, manganese metal, manganese dioxide, ferric ions, and many other substances such as yeast or even blood. A *catalyst* is a substance that, when added to a reaction mixture, participates in the reaction and speeds it up, but is not itself consumed in the reaction. The iodide ion is used as a catalyst in this demonstration.

Acknowledgment

Special thanks to Walter Rohr of Eastchester High School in Eastchester, NY for bringing this demonstration to our attention.

Materials for *Sudsy Kinetics* are available from Flinn Scientific, Inc.

| Catalog No. | Description |
|-------------|---|
| A0126 | Alconox, 4 lb |
| H0037 | Hydrogen Peroxide, 30%, 100 mL |
| H0028 | Hydrogen Peroxide, 6%, 500 mL |
| H0009 | Hydrogen Peroxide, 3%, 500 mL |
| S0084 | Sodium Iodide, 100 g |
| AP8599 | Hydrometer Cylinder |
| AP5429 | Demonstration Tray |
| AP4866 | Sudsy Kinetics—Chemical Demonstration Kit |

Consult your *Flinn Scientific Catalog/Reference Manual* for current prices.